

The materials and information provided via this brochure are distributed internationally, but not all the products or services found on this brochure are necessarily available in your country or area. Please contact local office or agents to make sure of the availability of the relevant products or services.

Mindray Medical Netherlands B.V.
Drs.W.van Royenstraat 8, 3871 AN Hoevelaken, The Netherlands
Voice: (31-33) 254-4911
Fax: (31-33) 253-4280

Mindray Medical Italy S.r.l.
Via Leonardo DA Vinci, 158, 20090 Trezzano Sul Naviglio (MI), Italy
Voice: (39-02) 5737-401
Fax: (39-02) 5523-2018

Mindray Medical España S.L.
Avenida Manoteras 38, Edificio B, Bajo B008. 28050-Madrid Spain
Voice: (34-91) 392-3754
Fax: (34-91) 392-3756

Mindray Medical Germany GmbH
Goebelstrasse 21, 3. OG, 64293 Darmstadt
Voice: (49 6151) 3910 0
Fax: (49 6151) 3910 300

Mindray (UK) Limited
Mindray House, Kingfisher Way, Hinchingsbrooke Business Park,
Huntingdon, Cambs, PE29 6FN
Voice: (44) 1480-416-840
Fax: (44) 1480-436-588

Mindray Medical France SARL
Europarc Creteil 1 allée des cerisiers 94000 CRETEIL, France
Voice: (33-1) 4513-9150
Fax: (33-1) 4513-9151

Mindray Medical Sweden AB
Rissneleden 136, 174 57 Sundbyberg, Sweden
Voice: (46) 8 555 54 121

Mindray Medical Rus Co.Ltd
2 Zvenigorodskaya Street, 13, bld 41, 123022 Moscow, Russia
Voice: (7-499) 553-6036
Fax: (7-499) 553-6039

Mindray Medical Colombia S.A.S
Room No.702, Torre UNIKA,Carrera 9 No.77-67, Bogota, Colombia
Voice: (57-1) 3130892 / 3210916

Mindray do Brasil – Comércio e Distribuição de Equipamentos Médicos Ltda.
Rua Tavares Bastos, 329, Perdizes, São Paulo - SP;
CEP 05012-020
Voice: (55-11) 3124-8026
Fax: (55-11) 3078-8035

Mindray Medical Mexico S. de R.L. de C.V.
Félix Parra # 175 Colonia. San José Insurgentes Sur Delegación Benito Juárez 03900 México, D. F.
Voice: (52-55) 5661-9450 / (52-55) 5662-6620
Fax: (52-55) 5662-6597

PT.Mindray Medical Indonesia
The East 11th Floor, Unit 08 Jl. Dr. Ide Anak Agung Gde Agung Kav. E 3.2 No.1 Jakarta 12950
Voice: (6221) 2902-7280
Fax: (6221) 2902-7283

Mindray Medical Egypt
Ground Floor -11 El Imam Metwally El Sharawi Street,
Sheraton Heliopolis - El Nozha, Cairo, Egypt
Voice: (2) 02 22691047
Fax: (2) 02 22691046

Mindray Medical Australia
ADD: 19 Harker St, Burwood VIC 3125, Australia
Voice: (61-3) 9830 6610
Fax: (61-2) 9550 2841

Mindray Medical Nigeria

Mindray Medical India Pvt. Ltd. (Delhi Office)
16th floor, Building No. 9-B, DLF Cyber City, DLF Phase- 3, Gurgaon- 122002, Haryana
Voice: (91-124) 4632-488
Fax: (91-124) 4632499

Mindray Medical India Pvt. Ltd. (Mumbai office)
B-404, City Point, Andheri Kurla Road, Andheri East, Mumbai, 400 059, India
Voice: (91-22) 4020-0000
Fax: (91-22) 4020-0011


Mindray Medikal Teknoloji Istanbul Ltd. Sti.
KÖYALTİ MEVKİİ 29 EKİM CAD.N.1 KUYUMÇUKENT 2.NO.PLZ.3.KT. D.1- 9 34520 BAĞÇELİEVLER/YENİBOSNA, İSTANBUL, TÜRKİYE
Voice: (90-212) 482 08 77
Fax: (90-212) 482 08 78

Mindray Medikal Teknoloji Istanbul Ltd. Sti. Branch Office Ankara
TURAN GÜNEŞ BULVARI ALEKSANDER DÜBÇEK CAD. NO:18/A ÇANKAYA/ANKARA, TÜRKİYE
Voice: (90-312) 441 96 92
Fax: (90-312) 441 96 93

Mindray Medical Vietnam Co.,Ltd.
Unit 2, 4th Floor, "B" Building, Broadway Office Park,
102 Nguyen Luong Bang Blvd., Tan Phu Ward, Dist. 7, HCMC, Vietnam
Voice: (84) 8 5413-6970
Fax: (84) 8 5413-6971

Mindray Medical Thailand Limited
1768 Thai summit Tower Room 407, 4th floor, New petchaburee Road, Kwaeng Bangkokpi, Khet Huai Khwang, Bangkok 10320, Thailand
Voice: (66) 2652-7288
Fax: (66) 2652-7266

Mindray Medical Korea
KNI Tower 9F 823-19 Yeoksam dong, Seoul-ganum, Korea
Tel: 008225688040
Fax: 008225688042

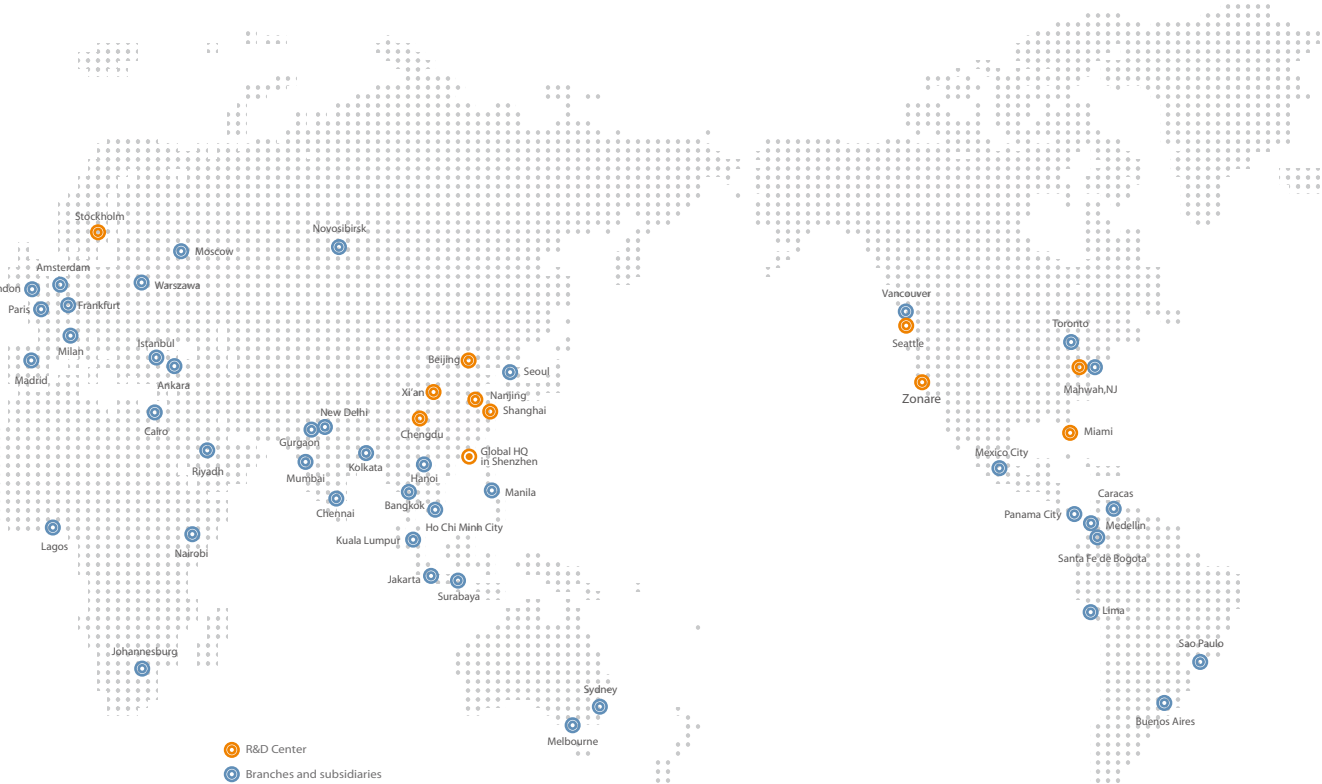

Mindray Building, Keji 12th Road South, High-tech Industrial Park, Nanshan, Shenzhen 518057, P.R. China Tel: +86 755 8188 8998 Fax: +86 755 26582680 E-mail: intl-market@mindray.com www.mindray.com	mindray healthcare within reach are registered trademarks or trademarks owned by Shenzhen Mindray Bio-medical Electronics Co., LTD. © 2015 Shenzhen Mindray Bio-Medical Electronics Co., Ltd. All rights reserved. Specifications subject to changes without prior notice. P/N:ENG-new SynoVent E-Series-210X285X8PX20160218
	mindray





new
SynoVent E-Series
Ventilator



mindray
healthcare within reach

Support for All Patients

- With an adjustable tidal volume down to 2 ml the SynoVent E-series ventilator is suitable for all patients categories from neonates to adults.
- For improved patient safety E-series ventilator utilises a proximal neonatal flow-sensor with minimal dead-space (<0.75ml) and low resistance (0.9cmH₂O @ 10L/min) to ensure total precision and control throughout the ventilation period.

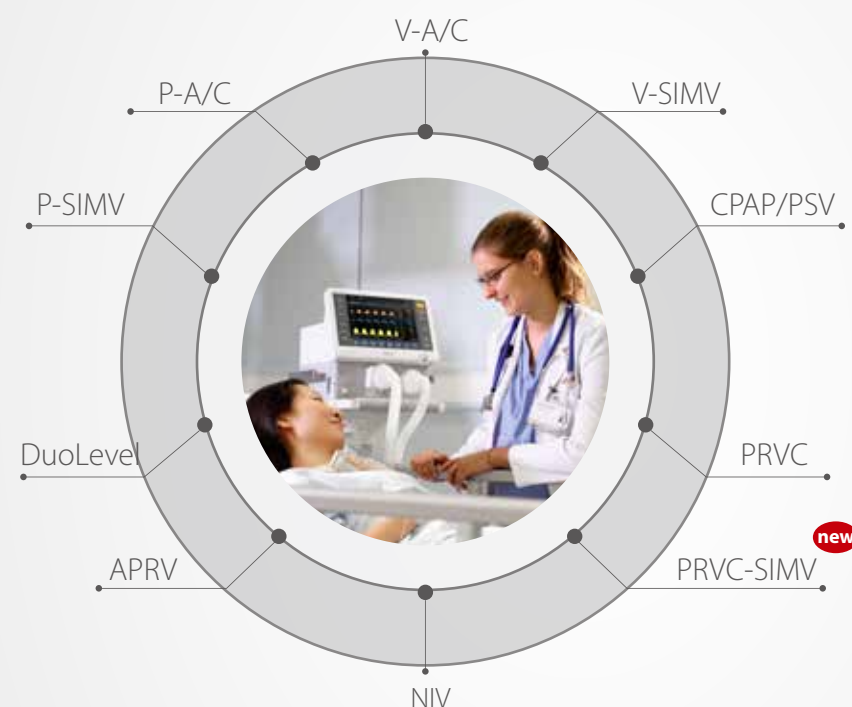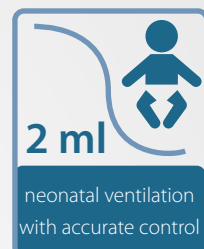

Enhanced Decision Making

- The continuous monitoring of up 31 parameters, including Work Of Breathing (WOB) and Intrinsic PEEP (PEEPi), ensures that all clinical decisions are evidence based and can be made with confidence.
- Both the side-stream and main-stream CO₂ modules are suitable for invasive and non-invasive ventilation. Furthermore the Volumetric CO₂ measurement capability reduces the requirement for repeated blood gas analysis and provides measurement of dead space and CO₂ elimination.
- For further support procedures can be initiated to measure Rapid Shallow Breathing Index (RSBI), Negative Inspiratory Force (NIF), and P0.1 for weaning process.




High-flow O₂ Therapy

- The E-series ventilator supports High-flow O₂ therapy ensuring delivery of precise oxygen concentrations at a user adjustable flow rate, optimal humidity and temperature.
- Delivery of High-flow O₂ therapy is shown to reduce CO₂ retention and improve hypoxia conditions which may improve prognosis.
- Switching between mechanical ventilation and O₂ therapy is achieved by enabling the O₂ therapy mode and changing the mask or nasal cannula of O₂ therapy, without the need to move the patient or change the ventilation equipment.

Network Compatibility

Using the Beneview patient monitor series the SynoVent E-series is able to connect seamlessly to both HIS and CIS systems with ease. And with an extensive number of communication ports patient data can be downloaded quickly for transfer and review.

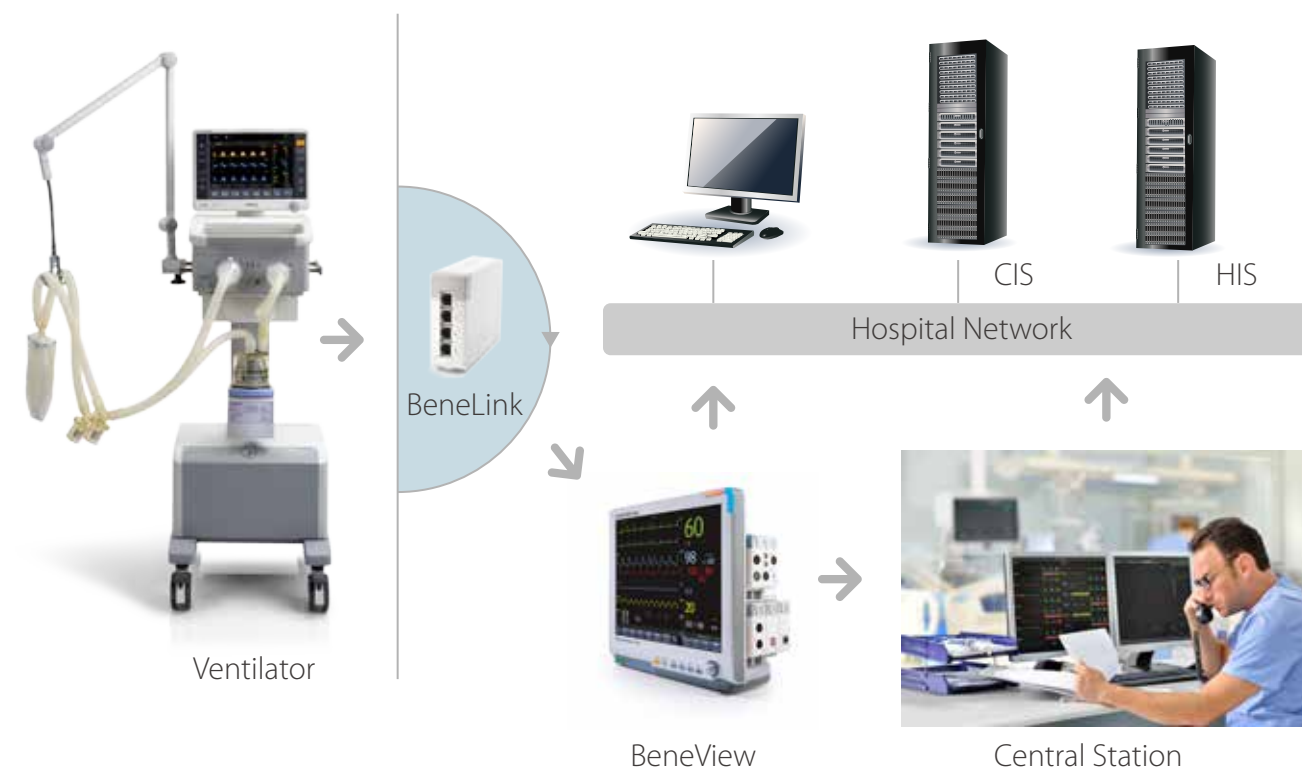



Widget: TVe/IBW

- For added patient safety the SynoVent E-series ventilator offers a simple and effective TV Setting prompt function.
- Each time changes are made to the set Tidal Volume value an on-screen prompt appears displaying the TVe / IBW. This feature may help in assisting the clinician to provide patients with improved, individualized, lung protective ventilation strategies. Doctors can even cross check it later on the trend page. This serves as a better guideline for patient lung protective strategies.



PV-tool

- PEEP is known to prevent atelectasis, improve gas exchange and reduce lung injury associated with repeated alveolar collapse / inflation.
- The setting of appropriate PEEP is one of the main goals for protective lung ventilation. With Mindray P-V Tool, the optimal PEEP for patients can be determined based on the characteristic points on the static P-V curve displayed by the SynoVent E-series ventilator .



Intuitive Display

With a large, user configurable, full colour touch screen the SynoVent E-series ventilator ensures that the data you want to view is clearly displayed at all times. Features include:

- Upto 4 waveforms display.
- Upto 3 loops display.
- Loop / waveform split screen display.
- Single page display of user selected parameters.
- Patients spontaneous breaths clearly highlighted in different colour for clear identification.



Intellicycle®

- Intellicycle® in E-series ventilator is designed to enhance patient / ventilator synchronization and improve patient comfort by automatically adjusting ventilatory support in accordance with changes in patients breathing pattern.
- The Intellicycle® feature is applicable to all patient categories and is useful in reducing the clinical workload of healthcare personnel.

Smart Apnea Backup

- Apnea backup ventilation will start to ensure patient safety by avoiding suffocation if the patient does not trigger any breath after the apnea time (user selectable).
- When patient able to trigger two consecutive breath within the apnea time, the ventilator will withdrawn back to the original mode.

Smart Suction

The SynoVent E-series ventilator is a supplied with a totally automatic suction protocol designed to simplify clinical workflow including:

- O₂ flush pre and post suction therapy.
- Alarm silence feature during suction technique.



Intelligent Standby for C3

For added safety the intelligent stand-by function of the C3 compressor ensures the patient is fully ventilated at all times by automatically maintaining ventilation even if the central gas supply fails.

